
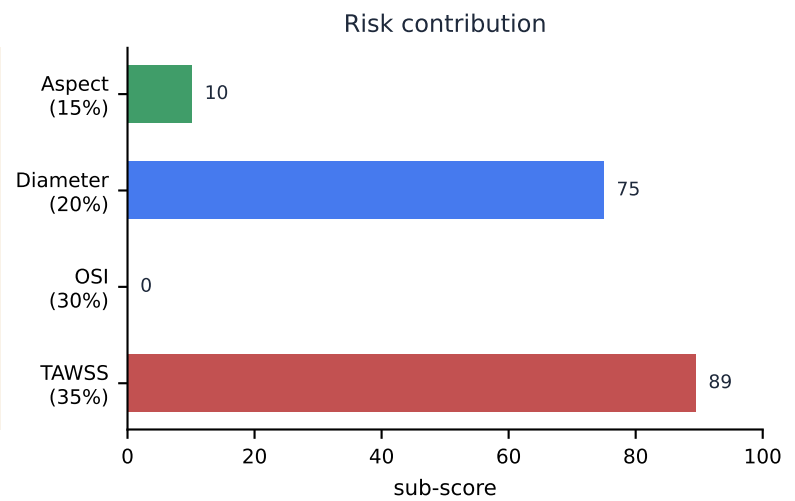


NeuroFlow CFD Analyst

Clinical Hemodynamic Assessment

Case PT-2026-0101 · generated 2026-08-05 11:27

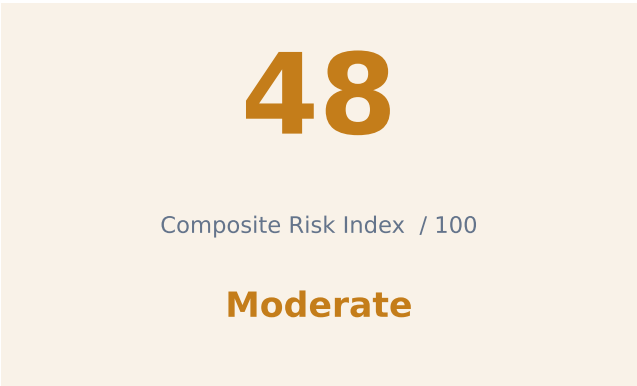

Hemodynamic parameters

Parameter	Value	Threshold	Status
TAWSS — aneurysm dome	0.292 Pa	< 0.4 Pa	FLAGGED
TAWSS — parent artery	3.291 Pa	reference	normal
OSI — aneurysm dome	0.008	> 0.2	normal
RRT — aneurysm dome	8.64 Pa ⁻¹	> 3.0	FLAGGED
ECAP — aneurysm dome	0.187 Pa ⁻¹	> 1.0	normal
NWSS (sac / parent)	0.089	—	normal
LSAR (relative, <10% parent)	82.8 %	—	normal
LSAR (absolute, <0.4 Pa)	86.0 %	—	normal

Morphology (measured from the reconstructed surface)

Max diameter	Neck diameter	Aspect ratio	Dome-to-neck
8.00 mm	6.90 mm	0.88	1.16
Volume	Surface area	Non-sphericity	
240.7 mm ³	174.1 mm ²	0.147	


Provenance

Solver: OpenFOAM ESI v2412 · Mesh: cells: 511015 · Time = 1.2e-05
WSS converted from OpenFOAM kinematic units (m²/s²) to Pascals by multiplying by rho = 1060.0. Validated against the analytic Poiseuille solution $\tau = 4 \cdot \mu \cdot Q / (\pi \cdot r^3)$.

NeuroFlow CFD Analyst

Zone Analysis and Model Explainability

Case PT-2026-0101

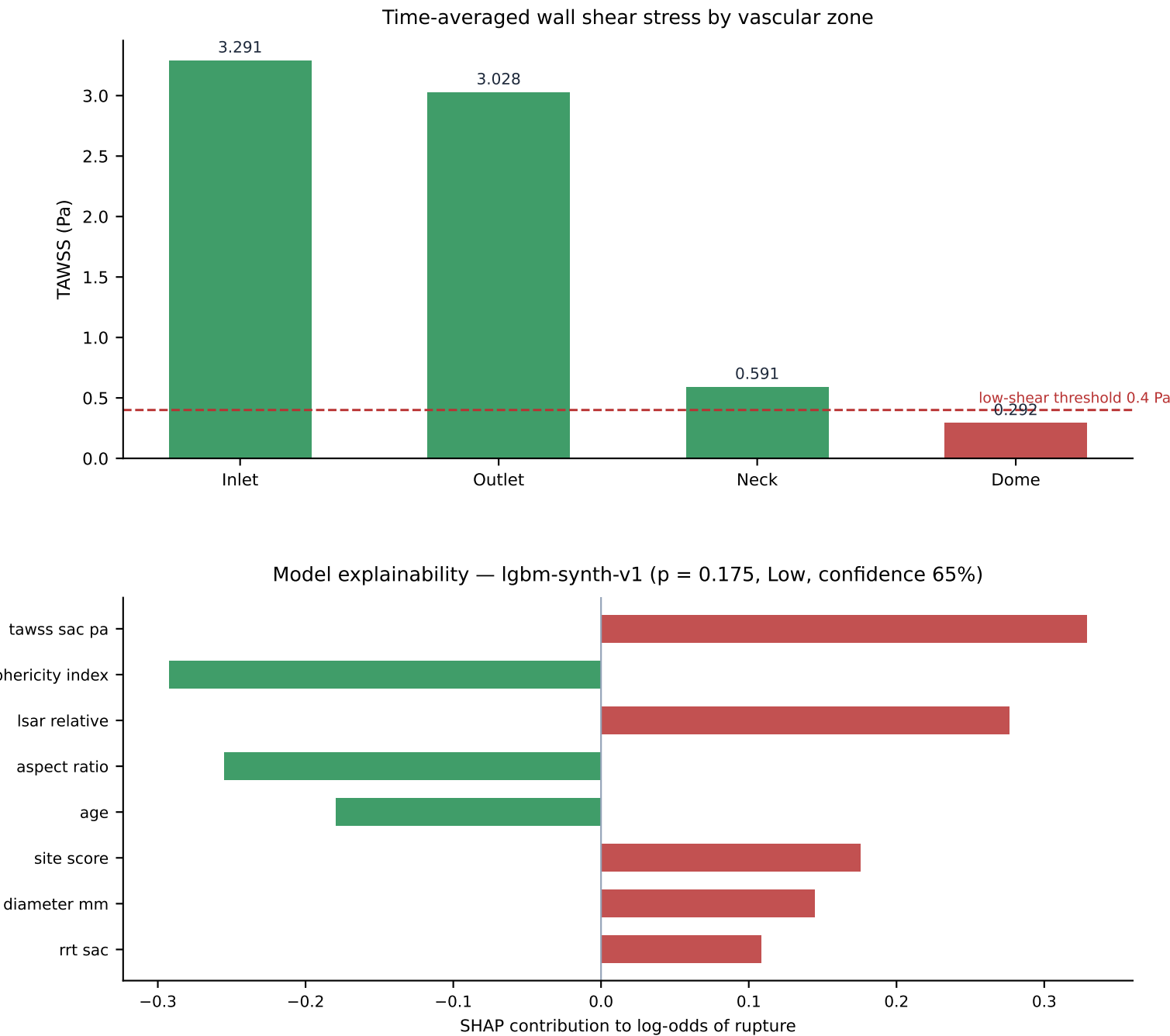

Assessment

Computational fluid dynamics analysis of case PT-2026-0101. Time-averaged wall shear stress at the aneurysm dome measured 0.292 Pa against 3.291 Pa in the parent artery (normalised WSS 0.089), with an oscillatory shear index of 0.008. 82.8% of the sac surface lies below 10% of the parent-artery shear (LSAR). Relative residence time at the dome is 8.64 Pa⁻¹. Sac morphology measured from the reconstructed surface: maximum diameter 8.00 mm, neck 6.90 mm, aspect ratio 0.88, volume 240.7 mm³. Composite Risk Index 48/100 (Moderate).